



KISTLER

measure. analyze. innovate.

Biomechanics

**Measuring Systems
for Performance
Diagnostics and Gait
and Balance Analysis
in Sports and Medicine**



Kistler – your partner for performance improvement and rehabilitation

Escalating athletic competitiveness can only be achieved through continuously improved training methods. An accurate knowledge of take-off forces at ski jumps or in athletics, and of force development involved in weightlifting or the martial arts helps set records time and again.

Our sense of balance is the product of a complex control system. Kistler force plates can be used for precise monitoring of shifts in weight and examination of the effects of diseases or medication.

Rehabilitation and prosthesis enhancement require very precise information on body movements, with the forces involved in walking and running playing a key role. Kistler force plates provide a precise and reliable record of these forces in every direction.

In addition to biomechanics measurement technology, as a Swiss company we also offer special sensors for measuring pressure,

force, torque and acceleration, as well as monitoring systems for mechanical production, the development and monitoring of internal combustion engines, automotive engineering and plastics processing.

Kistler's core competency is the development, production and use of sensors for measuring

- Pressure
- Force
- Torque
- Acceleration

With the aid of the company's expertise and electronic systems, measurement signals can be conditioned and used to analyze, control and optimize physical and other processes, and to boost product quality in manufacturing industry.

Year after year Kistler invests 10 % of its sales in R&D to facilitate technically innovative yet cost-effective state of the art

solutions. With a combined workforce of around 800, the Kistler Group is the market leader in dynamic measurement technology.

22 group companies and more than 30 distributors worldwide ensure close contact with the customer, individualized application engineering support and short lead times.

Kistler places a great deal of value on the exchange of know-how and close collaboration with leading research institutes, hospitals and sports performance centers worldwide.

This is the only way of developing reliable measuring instruments that fully meet the most stringent requirements.



Kistler is therefore a founder member of the International Society of Biomechanics (ISB), whose activities it supports as principal sponsor.

www.isbweb.org



Urs Kolly is a seven times Paralympic winner and has won many medals in European and World Championships

Kistler is a world champion in various measurement disciplines! An exact knowledge of force development helps achieve top-class performances time and again

Contents

Measurement technology in practice	Motion and gait analysis	4	
	Sport and performance diagnostics	6	
	Ergonomics and general biomechanics	8	
The piezoelectric effect		10	
Measuring chains with piezoelectric sensors		11	
Advantages of piezoelectric force plates		12	
Worthwhile investment without compromises		13	
Sensor portfolio biomechanics – general		14	
Product overview (applications)		16	
Product details/product range			
Measuring:	Force plates	17	
Connecting:	Cables/connectors	20	
Amplifying:	Charge amplifiers/control units	21	
Analysis:	Software	22	
Why Kistler force plates?		23	

KISTLER
measure. analyze. innovate.

measure.

Kistler develops and produces high quality measurement technology. Its core competence is sensor-related expertise.

analyze.

The analysis of sensor data calls for in-depth knowledge of the processes involved. Kistler has acquired this in major fields of application and is able to offer complete analysis systems.

innovate.

Kistler products are developed in close cooperation with key customers and universities. This inspires product innovation and technological progress.

Measurement technology in practice



Gait analysis laboratory with infrared cameras and force plates Type 9286AA



Special-sized force plates, e.g. for stairs (270x500x35 mm, top right) or for experiments with small animals (120x200x35 mm, top), are developed to customer requirements

Ever since Kistler force plates were first introduced back in 1969, they have consistently proved their worth as precise and reliable measuring instruments. They are a routinely used and indispensable tool in the motion analysis laboratories of leading specialists in biomechanics throughout the world. Their versatility and long service life have always made them a sound investment.

Their high sensitivity and extremely wide measuring range permitting one and the same plate to measure with equal precision forces both large and small. Their simple and dependable interfaces guarantee maximum long-term compatibility with other systems employed in motion analysis.

Fine resolution for gait analysis

Rehabilitation and the improvement of joint implants and prostheses demand a very precise understanding of body movements. The forces involved in walking and running play a key role here.

As well as measuring ground reaction forces precisely and reliably, Kistler force plates exploit inverse dynamics to provide a sound basis for calculating forces and moments. The high level of accuracy of these force plates detects minute changes and asymmetries in the gait. They therefore help practitioners reach correct diagnoses, adjust prostheses and document the rehabilitation process.

Motion and gait analysis





Source: Royal Veterinary College, University of London, UK

Maximum precision for equilibrium analysis

The human sense of balance is the product of a complex control system. In addition to precise force measurement, the special design principle of some Kistler force plates allows very accurate monitoring of shifts in weight and examination of the effects and progress of diseases or medication.

Kistler force plate Type 9286AA was developed specifically for gait and balance analysis. Even under a heavy load it permits measurement of extremely small changes in force and vibrations that can yield information about the condition of the neuromuscular motion control system.

Analysis of the gait of horses on six force plates Type 9287BA installed on the ground

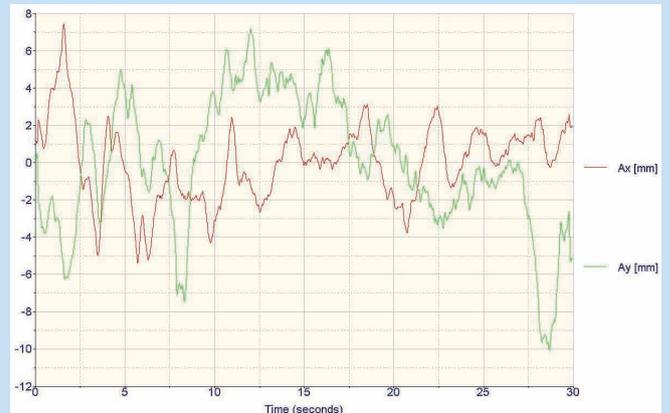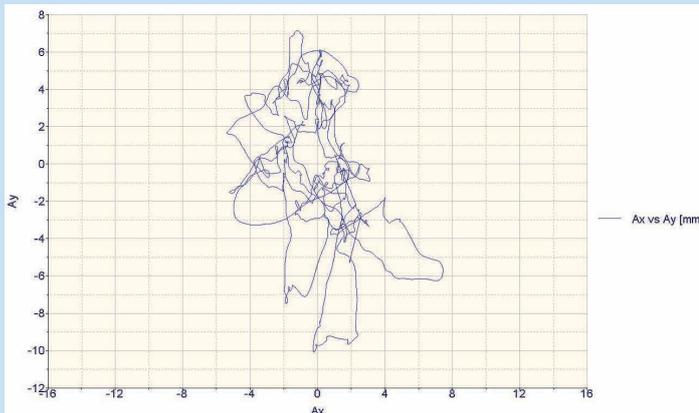

Equilibrium analysis through accurate determination of the center of pressure (COP); left: view of Ax against Ay; right: Ax and Ay against time



Force plate Type 9286AA is characterized by its very high COP accuracy, low overall height of just 35 mm and ease of installation on the floor

Measurement technology in practice

In sports, winning performances can only be achieved through continuously improved methods of highly controlled training. By analyzing the forces that athletes generate at various points, it is possible not only to draw conclusions relating to performance at a particular moment but also to record differences in the efficiency of various techniques and to modify the training program accordingly.

A precise knowledge of force development consistently helps prepare the way for record levels of performance. This applies equally well to technically demanding sports like golf, climbing and ski-jumping, but also to sports in which force (e.g. for explosive force) or energy (e.g. for endurance) need to be maximized or optimized.

Versatility means flexibility of application


The high preload and linearity of their piezoelectric sensors allow Kistler force plates to measure very accurately even under high initial loads. They can therefore be installed in any position and equipped with additional equipment or coverings without their precision or zero point being affected. The high natural frequency, extremely wide usable measuring range and high load capacity of the plates make it possible to capture even highly dynamic processes in the martial arts, Olympic lifting or athletics.



Force plate Type 9287BA for analysis of running and jumping in athletics



Ski jump with customized 3-component force plates for measuring forces on the jump-off platform

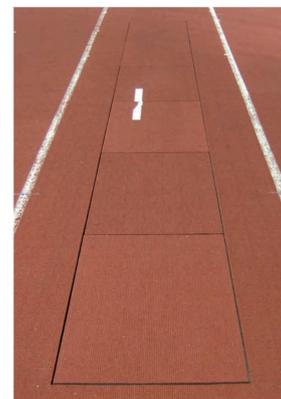

Sports and performance diagnostics

Measuring take-off power during natural movement

Kistler force plates provide an ideal instrument for objective measurement of take-off forces in performance diagnostics and biomechanics.

Quattro Jump is a complete mobile analysis system for testing the performance level of the legs in relation to power, coordination and stamina. In contrast to isokinetic systems or other jump tests the athlete moves completely naturally in a way that reflects his or her condition and coordination. Multiple jumps can be analyzed just as effectively as single jumps in order to provide objective, reliable and repeatable measurement results.

The reasonably priced easy-to-understand Quattro Jump quickly and accurately provides sports medicine physicians, trainers and athletes with the information needed to optimize training.

It is useful for:

- Performance and fitness testing and benchmarking
- Control and monitoring of training
- Matching the intensity of training to an injury



Source: Prof. Dr. Erich Müller, Interfaculty Department of Sports and Movement Science, University of Salzburg, Austria

Complete 3D motion analysis with determination of forces and moments in the knee joint during carved and conventional ski turns with the aid of four customized multicomponent force plates. The dynamometers are mounted between ski and binding/boot



Quattro Jump Type 9290AD consists of a portable force plate developed for measuring take-off force. For different jumps it measures the take-off force, which is then analyzed on the connected computer

Measurement technology in practice

Kistler provides a broad range of piezo-electric sensors for applications of every kind in research and industry. They are all characterized by an extremely wide measuring range, high linearity and stability and very compact size. This impressive selection is supported by a sales and service network providing vital expertise throughout the world.

Extensive portfolio of sensors for biomechanics

Our product line includes over 1 000 different force, acceleration, pressure and torque sensors for measuring tasks of all kinds. Force and movement are key factors in obtaining a clearer picture of the complex processes involved in biomechanics and a more precise definition of the properties of the materials. Because of Kistler's decades of accumulated experience in biomechanics, our experts are able to provide you with reliable and objective advice in selecting extremely versatile sensors and systems best suited for your special measurement task or application.



Medical treadmill ergometer



Instrumented horse treadmill with 18 Kistler force sensors

Instrumented treadmills

An instrumented horse treadmill was developed in the Sports Medicine Performance Center of the Equine Veterinary Clinic of the University of Zurich and constructed in conjunction with Kistler and the horse treadmill manufacturer Graber AG. The treadmill has 18 specially manufactured Kistler sensors and one 18-channel charge amplifier.

The "H/P/Cosmos Gaitway II" medical treadmill ergometer is based entirely on measurement technology and analysis software from Kistler and meets all of the legal standards relating to safety and medical devices. With force plates mounted in its bed, the treadmill allows quick and reliable clinical walking and galloping gait analysis.

High-sensitivity force sensors

Even such an extraordinary feat of measurement as determining the modulus of elasticity of a single trabecula of bone is possible with Kistler sensors. This structure between 1 and 2 mm in length and 50 ... 100 microns in thickness was deformed with the aid of a fine (50 micron) nylon thread in the natural network of trabecular bone. Force and deformation are measured to test the quality of the material forming the bone. High-sensitivity force sensor Type 9205 with a threshold of $<0,5 \cdot 10^{-3}$ N enabled measurement of the very low forces of up to 100 ... 1 500 mN.



Determination of modulus of elasticity of single trabecula of bone

Source: Prof. Dr. Edgar Stüssi, Institute of Biomechanics, ETH Zürich, Switzerland

Ergonomics and general biomechanics

Hand force measurement protects against strain

The forces we have to exert with our hands when working with heavy equipment and objects can be accurately measured over an extended period using a multicomponent hand force dynamometer. A knowledge of these forces for different activities helps in the diagnosis of occupational diseases and suggestion of preventive measures to reduce loads.

As with the use of force plates, the measured force vectors can be used to calculate the load on the spine and for other biomechanical evaluations.

The hand dynamometers were developed by the Institute for Occupational Safety and Health (BGIA) of the German employer's liability insurance association (Berufsgenossenschaften), using products including force sensors Type 9017B from Kistler.





Measurement of hand forces in different situations

Source: Berufsgenossenschaftliches Institut für Arbeitsschutz (BGIA), Germany



Determination of point of application of force and force vector on a climbing grip instrumented with two force sensors Type 9327A... during a National Climbing Championship in Singapore


Source: Prof. Franz Konstantin Fuss: Nanyang Technological University, Division of Bioengineering, Sports Engineering Research Team, Singapore; "Sports Technology", 2008

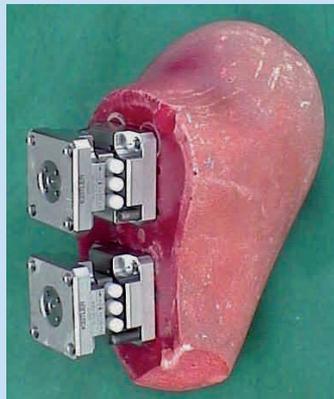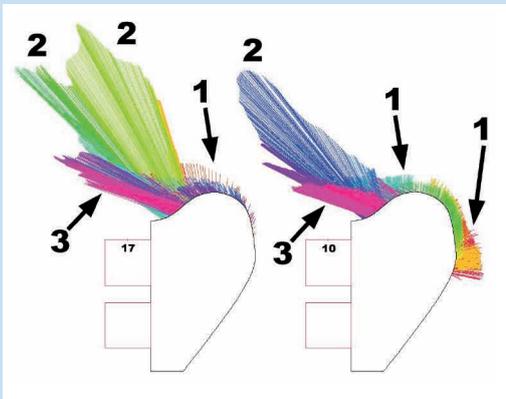

3-component force link 9327A for measuring tension and compression

The piezoelectric effect

Many crystals generate an electric charge when subjected to mechanical stress. This physical correlation has become known throughout the world as the piezoelectric effect. In 1950, Walter P. Kistler was granted the patent that heralded the breakthrough of wide industrial application for piezoelectric measurement. This method of measurement is the perfect answer to particularly extreme requirements in terms of geometry, temperature range and dynamics.

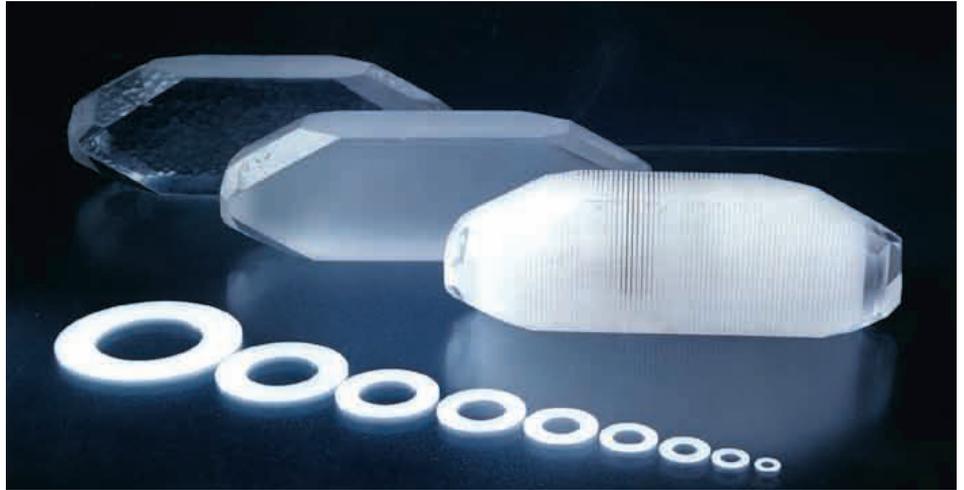

From quartz crystal to crystal washer

The piezoelectric effect – the prefix "piezo" comes from the Greek "piezein", to press – was discovered in 1880 by the Curie brothers. They found that the surfaces of certain crystals – including quartz – become electrically charged when the crystal is mechanically loaded. This electric charge is exactly proportional to the force acting on the crystal. It is measured in picocoulombs ($1 \text{ pC} = 10^{-12} \text{ coulombs}$).

Depending on the orientation of the polar axes of the crystal with respect to the applied force, two different effects relevant to biomechanics can be discerned

- Longitudinal
- Shear

Longitudinal effect

The charge produced by the longitudinal effect is developed on and can be collected from the surfaces to which the force is applied. Its magnitude Q in the case of the longitudinal effect depends only on the applied force F_x and not significantly on the dimensions of the crystal washers. The only way of increasing this charge is to connect several washers mechanically in series and electrically in parallel. The direction in which the crystal is sliced determines the properties and hence the application of the force link.


Piezoelements sliced to exhibit the longitudinal effect are sensitive to compression

forces and are therefore mainly suitable for simple, robust sensors for measuring forces.


Shear effect

As with the longitudinal effect, the piezoelectric sensitivity involved in the shear effect is independent of the size and shape of the piezoelement. The electric charge also develops on the loaded surfaces of the element in this case.


Shear-sensitive piezoelements are used for sensors measuring shear force, torque, strain and acceleration.

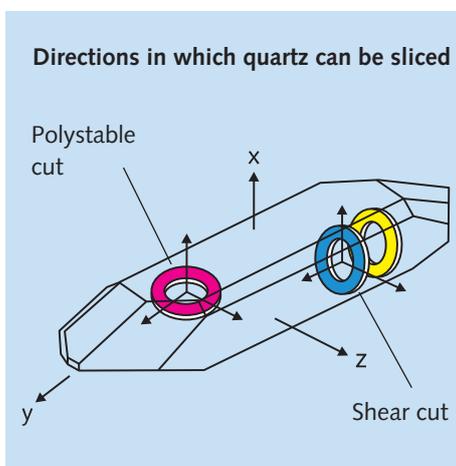

Directions in which quartz can be sliced

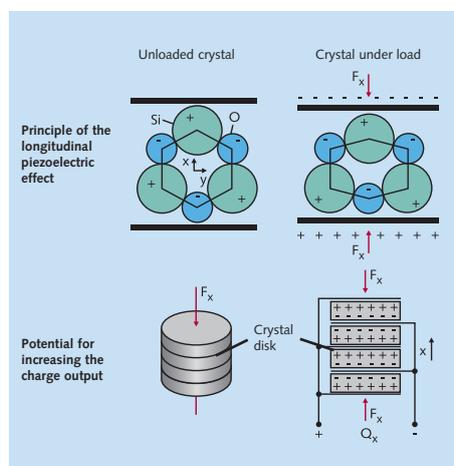

Principle of longitudinal effect

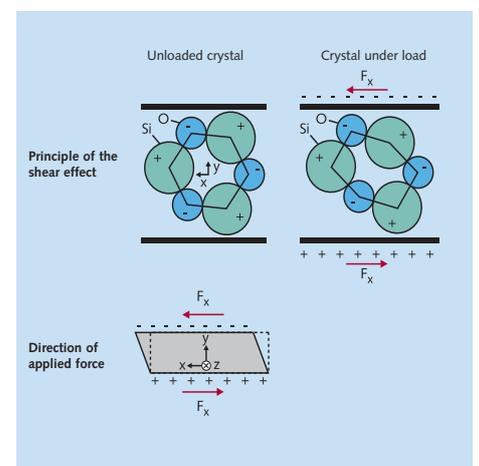

Principle of shear effect

Measuring chains with piezoelectric sensors

Types of sensor

Quartz washers with piezoelectric properties can be arranged in sensors to enable measurement of one or more force components or a torque vector. For use in biomechanics, Kistler offers the following sensors based on piezoelectrics:

- Single-component force sensors
- Multicomponent force sensors
- Multicomponent force plates and
- A wide range of other sensors for measuring strain, torque, pressure and acceleration


Single-component force sensors

Single-component force sensors, which are available in a wide variety of forms, are the design most suitable for measuring forces in a defined spatial direction.

Multicomponent force sensors

The piezoelectric measuring principle is also ideal for manufacturing multicomponent force sensors. The design of the sensor is similar to the single-component load washer. A pair of quartz washers sliced to exploit the longitudinal effect measures the normal component F_z , and one of each of a pair of washer sliced for the shear effect measures the two shear components F_x and F_y . As shear forces can only be transferred by friction, mounted multicomponent force sensors must always be under sufficiently high mechanical preload.

When mounted in dynamometers or force plates, multicomponent force sensors are usually employed in groups or three of four with the same sensitivity rather than individually. Clamping a load washer between

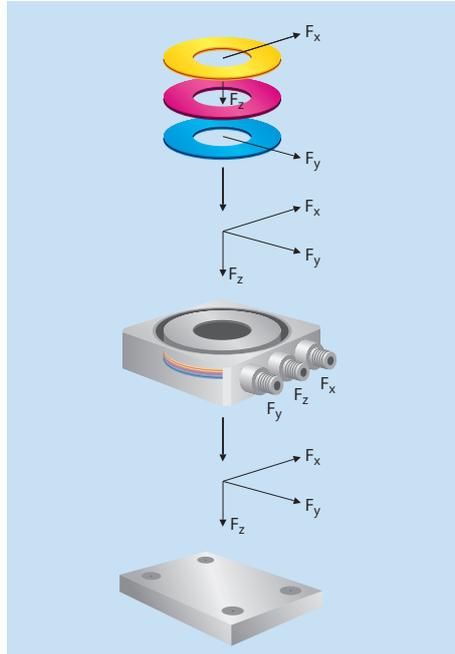

Principle of operation of multicomponent force plate

two special nuts gives rise to what is called a force link. This preloaded sensor can be used to measure tension and compression forces, for example in a linkage. Preloaded sensors are supplied calibrated and can be easily mounted ready to be used immediately for measurement.

Charge amplifiers

Charge amplifiers convert the charge output by a piezoelectric force sensor into a proportional voltage that can be used as an input variable for analysis or control systems. Most charge amplifiers from Kistler allow setting of the sensitivity and measur-

ing range of the sensor to allow use of one and the same sensor to measure very small through to very large quantities.

Calibration

Kistler sensors are calibrated for different measuring ranges prior to being dispatched from the factory. All of the relevant data is shown on the supplied calibration certificate. This ensures the output signal of the sensor can be converted accurately and reliably into the actual measurand (e.g. force). Kistler runs Swiss Calibration Service Laboratory No. 049, which is accredited to ISO 17025. The Kistler quality management system is certified to ISO 9001.

DAQ system and software

Kistler supplies a variety of powerful data acquisition systems with USB 2.0 ports or a PCI bus and its analysis software BioWare®, which is used to set charge amplifier parameters with great ease. All systems measure both highly dynamic processes and very small variables, and can also be used to acquire any analog signals. Connection cables and external control devices integrate force plates from Kistler into data acquisition and motion analysis systems from other manufacturers.

Configuration of typical measuring chain with Kistler DAQ System Type 5691A1

			
<p>Force plate with charge amplifier Type 9281CA</p>	<p>Connection cable Type 1759A...</p>	<p>DAQ-system (USB 2.0) Type 5691A1</p>	<p>Laptop (provided by user) with BioWare® software</p>

Advantages of piezoelectric force plates

Piezoelectric force, torque and strain sensors are very compact and rigid, and offer a measuring range of up to six decades, a high natural frequency and low interference sensitivity. They cover a wide temperature range, are overload-protected and offer long-term stability as well as freedom from fatigue. Piezoelectric sensors are ideal for almost all areas of application, particularly for the type of dynamic and highly sensitive processes encountered in biomechanics.

The operating principle of quartz crystal sensors in Kistler force plates means that compared with sensors with strain gages they offer decisive advantages, most of which are attributable to their comparatively high rigidity.

Sturdiness, overload protection and long-term stability

Piezoelectric force plates are very compact for an instrument with such a wide measuring range. Their rigidity makes them robust and gives them a high factor of safety against overload. Even after millions of load cycles and frequent thermal cycling the plates show no sign of fatigue or sensitivity drift. With correct use their service life is virtually unlimited.

High natural frequency and damping

The high rigidity of quartz crystal sensors is reflected in the very high overall stiffness of Kistler force plates, which results in a very high natural frequency and damping in all three directions of measurement. They are therefore particularly suitable for accurate measurement of highly dynamic processes such as the motion involved in sports.

High sensitivity over wide measuring range


Unlike strain gages, quartz crystal sensors have a wide measuring range of up to six decades, as their sensitivity, threshold and resolution are not linked to the stiffness of elastically deformable structures and hence to a limited measuring range. Irrespective of the measuring range and any initial load, piezoelectric sensors can also be used to measure very small forces, such in minute variations in a large base load.

Exact zero – No offset

With piezoelectric sensors the zero is re-determined prior to each new measurement process. Resetting physically tares the sensors by discharging any electric charge generated by factors such as an initial load. This automatically eliminates the effect of all static and slowly changing boundary conditions (such as mounting position, weight of supported elements and temperature fluctuations) without affecting the measurement accuracy.

Balance analysis and static measurements

Due to their low threshold, high sensitivity and very high accuracy Kistler force plates are ideal for balance analyses and other quasistatic measurements taken over a period of several minutes. However, piezoelectric sensors are not suitable for truly static measurements over several hours or days, as charge amplifiers always exhibit a

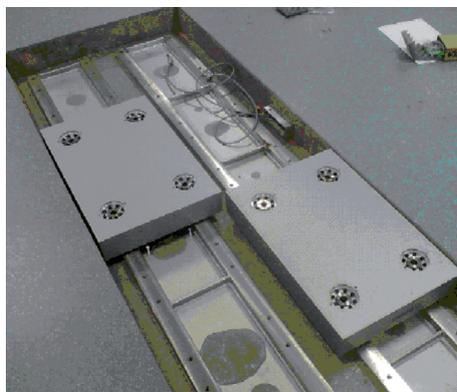

Easy installation and versatility: Two force plates Type 9281CA on frame with variable mounting positions

Benefits

Advantages of piezoelectric force plates over strain gage counterparts

- + Extremely wide measuring range of up to six decades
- + Sensitivity, threshold and resolution independent of measuring range
- + High rigidity and virtually displacement-free measurement
- + High natural frequency and damping
- + Overload protection, freedom from fatigue and long-term stability
- + Virtually unlimited service life
- + Low sensitivity to interference
- + Compact setup for wide measuring range

slight drift. Drift is defined as an undesirable change in output signal over an extended period of time that is not a function of the measurand. Charge amplifier drift is unavoidable because there is no material available offering infinite insulation.

Even the best transistors still exhibit minimal leakage currents of a few femtoamperes (10^{-15} A), and even the best insulators only have resistances of 10^{14} Ω . The maximum drift of Kistler charge amplifiers is approximately 0,03 pC/s, which in a force plate corresponds to a maximum unwanted change in the signal of $\leq \pm 10$ mN/s for F_z (or $\leq \pm 5$ mN/s for F_x, y).

Worthwhile investment without compromises

Kistler force plates measure extremely accurately and provide reliable, reproducible results to guaranteed specifications over their entire service life. They quickly pay for themselves, as they are built to last and cover a wide range of applications and extensive measuring range without compromising on measurement accuracy.

Easy installation and mobility


In contrast to conventional designs the force plate Type 9286A... for gait and balance analysis does not have to be installed on a frame; it can just be used on any flat surface without incurring additional costs. The plate's low overall height of just 35 mm

and weight of under 18 kg allows versatile, mobile application. For very dynamic applications Kistler force plates are mounted on an installation frame that is molded into the foundation to guarantee their full potential can be exploited. A lightweight aluminum sandwich cover plate keeps the force plates Type 9281C... and 9287B... relatively light at 16 or 25 kg. This makes installation and transport easier, so that Kistler force plates can be used not only in different mounting positions but also in various measurement locations.

Cost effectiveness through versatility

While force plates with strain gage sensors are limited by their elastically deformable structures to a particular, optimal measur-

ing range, Kistler force plates with quartz crystal sensors can accurately measure both very small and very large forces within their structural load capacity. Very small forces can also be measured under a high initial load. It is therefore not necessary to purchase different Kistler force plates to cover several measuring ranges. Moreover, quartz crystal sensors are characterized by extreme robustness and durability. Due to their high natural frequency in all three directions of measurement (f_{nx} , f_{ny} and f_{nz}), high natural damping and rapid response, Kistler force plates always represent even highly dynamic processes correctly.

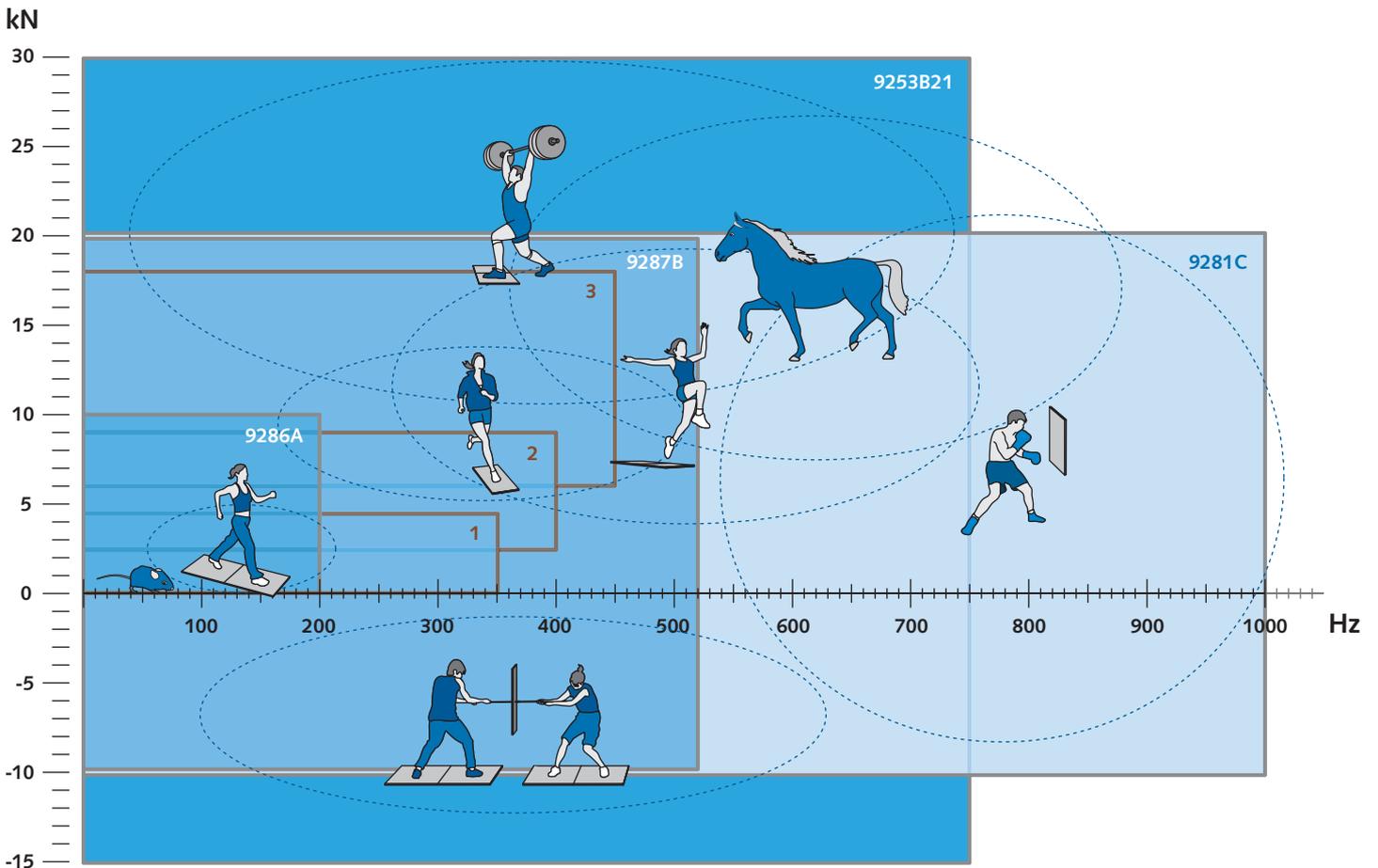

Kistler force plates are versatile and suitable – depending on maximum load capacity and natural frequency – for a wide variety of applications. In comparison three typical force plates with strain gage (4,5, 9 and 18 kN) and only limited application fields

y-axis: maximum load capacity
x-axis: lowest natural frequency from f_{nx} , f_{ny} and f_{nz}



Typical force plates with strain gages:


- 1 strain gage F_z max. 4,5 kN
- 2 strain gage F_z max. 9 kN
- 3 strain gage F_z max. 18 kN

Sensor portfolio for biomechanics – general

Force and movement are key factors in gaining a clearer idea of the complex processes of biomechanics and a more precise definition of material characteristics.

Kistler's range includes over 1 000 different force, acceleration, pressure and torque sensors for measuring tasks of all kinds.

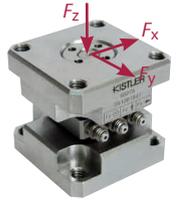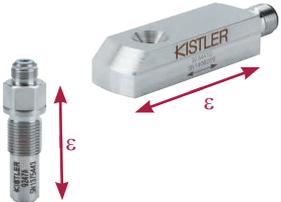Description		Range	Product
Force	<p>Force sensors Quartz sensors from Kistler can be used for direct and indirect measurement of forces in one or more directions.</p>	1 mN ... 20 MN	
	<p>Sensors for very small forces The special design of these sensors provides sensitivity 30 times higher than load washers.</p>	-20 ... 500 N	
	<p>Load washers Kistler robust load washers are extraordinarily versatile in application.</p>	0 ... 1 200 kN	
	<p>Miniature sensors The extremely small size, high natural frequency and integral connecting cable enable use under critical mounting conditions.</p>	0 ... 2,5 kN	
	<p>1-component force link The easily mounted preloaded force link is calibrated ready for immediate measurement of tension and compression forces.</p>	-20 ... 120 kN	
	<p>Shear elements These very flat and small elements allow extremely flexible mounting.</p>	-4 ... 4 kN	
	<p>Strain sensors (longitudinal/transverse) This type of sensor measures the strain of the structure in which it is mounted and provides indirect measurement of very high forces.</p>	-600 ... 600 $\mu\epsilon$, -1 500 ... 1 500 $\mu\epsilon$	

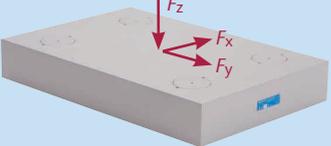	Description	Range	Product
Force	<p>Load washers Kistler 3-component load washers measure the three orthogonal force components independently.</p>	-2 ... 2 kN	
	<p>3-component force link The easily mounted preloaded force link is calibrated ready for immediate measurement of the three force components.</p>	-2 ... 2 kN	
	<p>With built-in electronics 3-component load washers with built-in charge amplifier and output voltage (5 V).</p>	-5 ... 5 kN	
	<p>Multicomponent force plate This plate uses 4 multicomponent load washers to enable measurement of the orthogonal forces, the moments and the point of application of the force.</p>	-10 ... 30 kN	
	<p>Kit Ready-to-connect kit for assembling multicomponent force plate.</p>	-20 ... 40 kN	
	<p>Torque and force Reaction torque sensor capable of measuring additional acting forces.</p>	F_z -10 ... 10 kN, -20 ... 20 kN	
Pressure	<p>Pressure sensors Kistler supplies technology for static and dynamic pressure measurement characterized by reliability, precision and flexibility.</p>	0,1 ... 10 000 bar	
Acceleration	<p>Accelerometers Kistler offers an extensive selection of different accelerometers covering extreme measurement requirements.</p>	3 μ g ... 100 000 g	

More detailed information and datasheets are to be found in the Product Finder on Kistler's website.

Product overview (applications)

Force plate						
Application	Type 9285	Type 9286A...	Type 9281C...	Type 9287B...	Type 9253B...	Type 9290AD
Take-off force in performance diagnostics	—	○	+	+	+	++
Mobile force measurement (field trials)	—	++	+	○	—	+
Highly dynamic force measurement in competition	—	○	++	++	+	—
Force measurement under difficult conditions (high load, point load, moisture, etc.)	—	+	+	+	++	—
Force measurement under very difficult conditions (highly dynamic, large/small forces, wide measuring range, etc.)	—	+	++	++	++	—
Gait and running analysis (varied and dynamic)	+	++	++	+	+	—
Gait and balance analysis with very accurate determination of COP, fixed installation	++	++	+	+	+	—
Clinical gait and balance analysis with very accurate determination of COP, portable, low installation costs	○	++	○	—	—	—

Sports
Research
Clinic


Key: Ideal ++ Very suitable + Possible ○ Unsuitable —

Necessary requirements	Fulfilled					
Very dynamic processes, wide useful frequency range			•	•	•	
Very high forces (F >10 kN)			•	•	•	
Very low forces (F <100 N), high sensitivity	•	•	•	•	•	•
Very accurate determination of COP	•	•				
Very heavy point loading of cover plate		•			•	
Waterproof (IP67)					•	
Glass cover plate	•					
Integral charge amplifier available		•	•	•		•

Measuring

Force plates

Multicomponent force plate

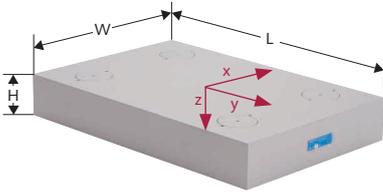

Type 9281C...

Technical data			Type 9281C...
Measuring range	F_x, F_y	kN	-10 ... 10
	F_z	kN	-10 ... 20
Overload	F_x, F_y	kN	-15/15
	F_z	kN	-15/25
Linearity	%FSO		< $\pm 0,5$
Hysteresis	%FSO		<0,5
Crosstalk	$F_x \leftrightarrow F_y$	%	< $\pm 1,5$
	$F_x, F_y \rightarrow F_z$	%	< $\pm 1,5$
	$F_z \rightarrow F_x, F_y$	%	< $\pm 1,0$ (inside sensor rectangle)
Rigidity	x-axis ($a_y = 0$)	N/ μm	≈ 250
	y-axis ($a_x = 0$)	N/ μm	≈ 400
	z-axis ($a_x = a_y = 0$)	N/ μm	≈ 30
Natural frequency	$f_n(x, y)$	Hz	$\approx 1\,000$
	$f_n(z)$	Hz	$\approx 1\,000$
Operating temperature range		$^{\circ}\text{C}$	0 ... 60
Weight		kg	16
L		mm	600
W		mm	400
H		mm	100
Degree of protection	EN 60529:1992		IP65

Properties

Extremely wide measuring range, excellent measuring accuracy, high natural frequency, versatile, threshold $F_z < 250$ mN.

Areas of application


This force plate is designed specifically for use in basic research, sports and gait analysis. The plate can be mounted in any position.

Accessories

DAQ system Type 5691A1
see data sheet for further details

Data sheet 9281C_000-156

Large multicomponent force plate

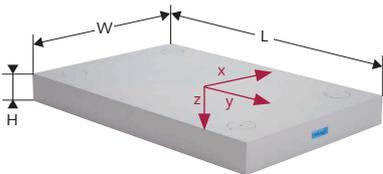

Type 9287B...

Technical data			Type 9287B...
Measuring range	F_x, F_y	kN	-10 ... 10
	F_z	kN	-10 ... 20
Overload	F_x, F_y	kN	-13/13
	F_z	kN	-13/25
Linearity	%FSO		< $\pm 0,5$
Hysteresis	%FSO		<0,5
Crosstalk	$F_x \leftrightarrow F_y$	%	< $\pm 1,5$
	$F_x, F_y \rightarrow F_z$	%	< $\pm 1,5$
	$F_z \rightarrow F_x, F_y$	%	< $\pm 1,0$ (inside sensor rectangle)
Rigidity	x-axis ($a_y = 0$)	N/ μm	≈ 150
	y-axis ($a_x = 0$)	N/ μm	≈ 200
	z-axis ($a_x = a_y = 0$)	N/ μm	≈ 30
Natural frequency	$f_n(x, y)$	Hz	≈ 750
	$f_n(z)$	Hz	≈ 520
Operating temperature range		$^{\circ}\text{C}$	0 ... 60
Weight		kg	25
L		mm	900
W		mm	600
H		mm	100
Degree of protection	EN 60529:1992		IP65

Properties

Extremely wide measuring range, excellent measuring accuracy, high natural frequency, versatile, threshold $F_z < 250$ mN, large dimensions.

Areas of application

This force plate is designed specifically for use in basic research, sports and gait analysis. The plate can be mounted in any position.

Accessories


DAQ system Type 5691A1
see data sheet for further details

Data sheet 9287B_000-159

Measuring

Force plates

Mobile multicomponent force plate

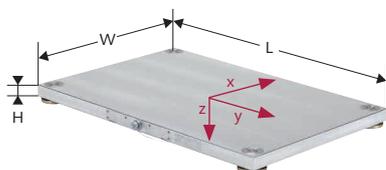

Type 9286A...

Technical data			Type 9286A...
Measuring range	F_x, F_y	kN	-2,5 ... 2,5
	F_z	kN	0 ... 10
Overload	F_x, F_y	kN	-3/3
	F_z	kN	0/12
Linearity	%FSO		< \pm 0,5
Hysteresis	%FSO		<0,5
Crosstalk	$F_x \leftrightarrow F_y$	%	< \pm 1,5
	$F_x, F_y \rightarrow F_z$	%	< \pm 2,0
	$F_z \rightarrow F_x, F_y$	%	< \pm 0,5 (inside sensor rectangle)
Rigidity	x-axis ($a_y = 0$)	N/ μ m	\approx 12
	y-axis ($a_x = 0$)	N/ μ m	\approx 12
	z-axis ($a_x = a_y = 0$)	N/ μ m	\approx 8
Natural frequency	$f_n(x, y)$	Hz	\approx 350
	$f_n(z)$	Hz	\approx 200
Operating temperature range		$^{\circ}$ C	0 ... 60
Weight		kg	17,5
L		mm	600
W		mm	400
H		mm	35
Degree of protection	EN 60529:1992		IP63

Properties

Excellent accuracy of center of pressure (COP), very wide measuring range, easy mounting, flexible, portable usage, threshold $F_z < 250$ mN.

Areas of application


This force plate is designed specifically for use in gait and balance analysis. Flexible, portable usage.

Accessories

DAQ system Type 5691A1
Walkway Type 9401B...
see data sheet for further details

Data sheet 9286A_000-158

Multicomponent force plate with glass top plate

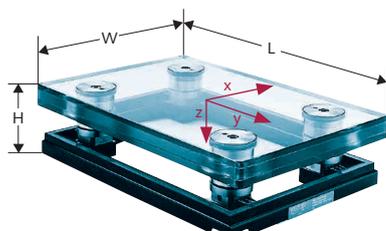

Type 9285

Technical data			Type 9285
Measuring range	F_x, F_y	kN	-2,5 ... 2,5
	F_z	kN	0 ... 5
Overload	F_x, F_y	kN	-3,75/3,75
	F_z	kN	0/12
Linearity	%FSO		< \pm 0,5
Hysteresis	%FSO		<1
Crosstalk	$F_x \leftrightarrow F_y$	%	< \pm 2
	$F_x, F_y \rightarrow F_z$	%	< \pm 2
	$F_z \rightarrow F_x, F_y$	%	< \pm 1
Rigidity	x-axis ($a_y = 0$)	N/ μ m	\approx 120
	y-axis ($a_x = 0$)	N/ μ m	\approx 115
	z-axis ($a_x = a_y = 0$)	N/ μ m	\approx 25
Natural frequency	$f_n(x, y)$	Hz	\approx 300
	$f_n(z)$	Hz	\approx 500
Operating temperature range		$^{\circ}$ C	-20 ... 50
Weight		kg	45
L		mm	600
W		mm	400
H		mm	150
Degree of protection	EN 60529:1992		IP65

Properties

Glass top plate allows recording of contact surface, wide measuring range, excellent measuring accuracy, excellent accuracy of center of pressure (COP), threshold $F_z < 10$ mN.

Areas of application

This force plate is designed for special gait and balance analysis applications. The glass plate allows simultaneous force measurement and photographic or cinematographic recording of the contact surface from below.

Accessories


Cable Type 1685B...
Charge amplifier Type 9865E...
see data sheet for further details

Data sheet 9285_000-157

Measuring

Force plates

Multicomponent force plate

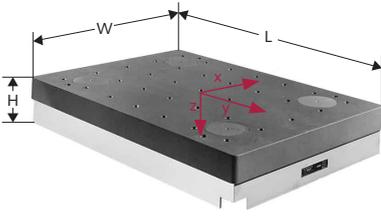

Type 9253B...

Technical data			Type 9253B11/12	Type 9253B21/22	Type 9253B23
Measuring range	F_x, F_y	kN	-10 ... 10	-15 ... 15	-12 ... 12
	F_z	kN	-10 ... 20	-15 ... 30	-12 ... 25
Overload	F_x, F_y	kN	-15/15	-20/20	-15/15
	F_z	kN	-15/30	-20/40	-15/30
Linearity	%FSO		<±0,5	<±0,5	<±0,5
Hysteresis	%FSO		<0,5	<0,5	<0,5
Crosstalk	$F_x \leftrightarrow F_y$	%	<±2	<±2	<±2
	$F_x, F_y \rightarrow F_z$	%	<±2	<±2	<±2
	$F_z \rightarrow F_x, F_y$	%	<±2	<±2	<±2
Rigidity	x-axis ($a_y = 0$)	N/μm	≈625	≈750	≈850
	y-axis ($a_x = 0$)	N/μm	≈650	≈850	≈750
	z-axis ($a_x = a_y = 0$)	N/μm	≈250	≈450	≈250
Natural frequency	$f_n(x, y)$	Hz	≈800, ≈750	≈580, ≈550	≈610, ≈570
	$f_n(z)$	Hz	≈850	≈720	≈570
Operating temperature range		°C	-20 ... 70	-20 ... 70	-20 ... 70
Weight		kg	40	90	85
L		mm	600	600	600
W		mm	400	400	400
H		mm	100	100	100
Degree of protection	EN 60529:1992		IP67	IP67	IP67

Properties

Heavy-duty aluminum or steel cover plate, waterproof, very wide measuring range.

Areas of application

Special applications under more difficult conditions: immersed in water, with heavy loads, problematic ambient conditions etc.

Accessories

Connection cable Type 1677AQ02sp
Charge amplifier Type 9865E...
see data sheet for further details

Data sheet 000-146

Quattro Jump, portable force plate system



Type 9290AD

Technical data			Type 9290AD
Measuring range	F_z	kN	0 ... 10
Overload	F_z	kN	15
Linearity	%FSO		<±0,5
Hysteresis	%FSO		<1
Natural frequency	$f_n(z)$	Hz	≈150
Operating temperature range		°C	0 ... 50
Weight		kg	21,6
L		mm	920
W		mm	920
H		mm	125
Sampling rate		Hz	500
Resolution	range 1	N/bit	1
	range 2	N/bit	9,2
Interface to the computer	Connector type		RS-232C
	Baudrate	kBaud	D-Sub 9f
			19,2 ... 115
Power supply	Battery, LRG, AMS	V	12x1,5
	External power supply	V	12
	Connector type jack-socket		D5,5/2,1
Battery lifetime (Alkaline) h		≈15	
Degree of protection	EN 60529:1992		IP50

Properties

User-friendly, portable usage, integral charge amplifier and digitization

Areas of application

Objective measurement of jump force, jump height and jump power, immediate feedback to optimize the training program.

Accessories

Battery (12x) Type 5.310.002
see data sheet for further details

Data sheet 9290A_000-160

Connecting

Connection cables

Cables




Technical data		Type 1681B...
Connection		Fischer 9 pin pos. – Fischer 9 pin pos.
Length	m	5 / 10 / 20
Diameter	mm	12,3
Degree of protection (EN 60529)		IP65




Technical data		Type 1685B...
Connection		Fischer 9 pin pos. – Fischer 9 pin pos.
Length	m	5 / 10 / sp (max. 20)
Diameter	mm	5,6
Degree of protection (EN 60529)		IP65




Technical data		Type 1686A...
Connection		Fischer Winkel 9 pin pos. – Fischer 9 pin pos.
Length	m	5 / 10 / sp (max. 20)
Diameter	mm	5,6
Degree of protection (EN 60529)		IP65

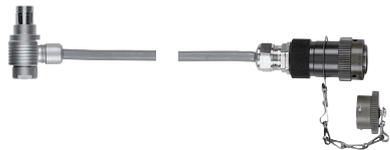

Technical data		Type 1757A...
Connection		Fischer angle 19 pin pos. – MIL 19 pos.
Length	m	10 / sp (max. 30)
Diameter	mm	5,6
Degree of protection (EN 60529)		IP63

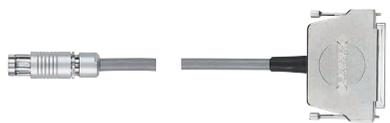

Technical data		Type 1758A...
Connection		Fischer 19 pin pos. – D-Sub37 neg.
Length	m	10 / sp (max. 30)
Diameter	mm	8
Degree of protection (EN 60529)		IP63

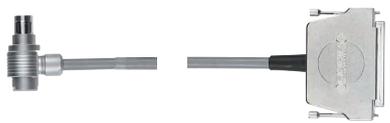

Technical data		Type 1759A...
Connection		Fischer angle 19 pin pos. – D-Sub37 neg.
Length	m	10 / sp (max. 30)
Diameter	mm	8
Degree of protection (EN 60529)		IP63



Technical data		Type 1760A...
Connection		Fischer 19 pin pos. – MIL 19 pos.
Length	m	10 / sp (max. 30)
Diameter	mm	8
Degree of protection (EN 60529)		IP63

Amplifying

Amplifier

8-channel charge amplifier



Technical data			Type 9865E...
Number of measuring channels			8
Measuring ranges	range 1	pC	±1 000
	range 2	pC	±5 000
	range 3	pC	±10 000
	range 4	pC	±50 000
Output voltage (selectable)		V	±5 / ±10
Output current		mA	<±5
Output impedance		Ω	10
Upper cut-off frequency		kHz	≈10
Lower cut-off frequency			
ranges 1 000 / 5 000 pC		s	≈10
ranges 10 000 / 50 000 pC		s	≈100
Error (all channels)		%FS	<1
Noise (at output)		mV _{rms}	<2
Drift		pC/s	<±0,07
Connections			
Measuring input			Fischer socket, 9 pin
Outputs (analog)			D-Sub 15 pin, neg.
Inputs (digital, TTL-S)			D-Sub, 9 pin, neg.
Power line (protection class I)			
Power line voltage (switchable)		V AC	230 / 115
Tolerance		%	-22 / 15
Frequency		Hz	48 ... 62
Power consumption		VA	≈25
Operating temperature range		°C	0 ... 50
Dimensions (WxHxD)		mm	236x151x225
Weight		kg	≈4
Degree of protection			
Electro-medical equipment			EN 60601-1, EN 60601-1-2
Safety requirements			EN 61010-1
EMC emission			EN 50081-1
EMC immunity			EN 50082-1

Properties

8-channel charge amplifier with 4 measuring ranges, remote controlled.

Areas of application

Universal, economically priced multichannel charge amplifier for Kistler force plates.

Accessories

Connection cable Type 1677A...
 Connection cable Type 1681B...
 Connection cable Type 1685B...
 Connection cable Type 1686A...
 BioWare® Type 2812A...
 see data sheet for further details

Data sheet 9865E_000-287

Control unit



Technical data			Type 5233A2
Number of measuring channels			8
Output signal		V	-10 ... 10
Connections			
Connection to force plate			MIL 1419 (19 pin)
Connection to PC (remote control)			D-Sub 37 neg.
Analog outputs			BNC neg.
Power (2 phases + ground)			IEC 320 C14
Power line, switchable		V AC	230 / 110
Power line frequency		Hz	48 ... 60
Dimensions		mm	170x126x55

Properties

External control unit for measuring components with built-in charge amplifier.

Areas of application

Manual control of Kistler multicomponent force plates with analog signal output.


Accessories

Connection cable Type 1757A10
 Connection cable Type 1760A10
 see data sheet for further details

Data sheet 5233A_000-150

Software

DAQ system for BioWare® 4.0

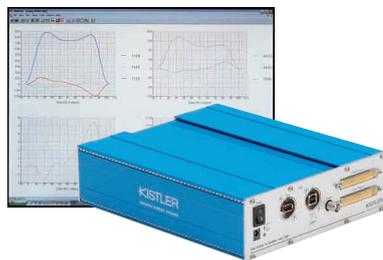

Technical data		Type 5691A1	
Dimensions	mm	208x65x250	
Weight (total)	kg	2,05	
Operating temperature range	°C	0 ... 50	
Power supply voltage			
Power supply	VDC	11 ... 15	
Power consumption	VA	6	
AD converter			
Number of channels		16	
Resolution (per channel)	Bit	16	
Input voltage range	V	±1, ±2, ±5, ±10 (software selectable)	
Sampling rate	S/s	0,6 ... 50 000 (software selectable)	
1 Force plate	max. @ 2 channels	kS/s	50
2 Force plates	max. @ 8 channels	kS/s	17
	max. @ 16 channels	kS/s	9,5
Connections USB 2.0			
USB In (uplink, to the PC)		USB Type B, female	
USB Out (downlink, free)		USB Type A, female	
Force plate 1/2			
Input voltage range (max.)	V	D-Sub37, male ±15	
External trigger	Type	BNC neg.	
Input voltage pull-up resistance			
10 kΩ on ±5 V			
max.	VDC	12	
high or input open	VDC	>3,6	
low	VDC	<0,6	
Trigger modus	standard	rising edge	
	software selectable	falling edge	

Properties

USB 2.0 ensures easy installation, remote control with integral charge amplifier, Powerful data acquisition and signal processing system, versatile data analysis and filter, external and internal trigger. Includes BioWare® 4.0 software.

Areas of application


Suitable for measuring both slow and highly dynamic processes as well as very small measurands. Suitable for basic research, sports science, gait analysis, ergonomics, etc.

Accessories

Connection cable Type 1758A...
Connection cable Type 1759A...
BioWare® Type 2812A...
see data sheet for further details

Data sheet 5691A_000-633

BioWare® 4.0

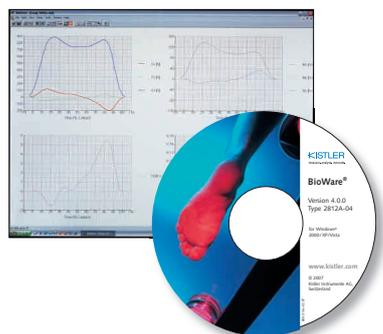

Properties

Very versatile, easy to use Windows® software specially designed to work with Kistler force plates in the various fields of biomechanics:

- Biomechanics research
- Gait analysis (rehabilitation, orthopaedics, prosthetics)
- Sports (jump force, impact, training)
- Neurology (posturography, balance, microvibrations)
- Ergonomics, industry (shoe development, material testing, safety, loading)

BioWare contains all data acquisition, signal conditioning and analysis of force plate data.

Areas of application

Suitable for use in basic research, gait analysis, sports, neurology, ergonomics and industry.

Accessories

BioWare® data Type 2812A...
acquisition system
see data sheet for further details

Data sheet 2812A_000-370

Windows® is a registered trademark of Microsoft Corporation.

Why Kistler force plates?



The unique piezoelectric measuring system offers numerous advantages over force plates with strain gages, and makes Kistler force plates a cost-effective investment for biomechanical and medical laboratories.

Advantages of Kistler force plates

- Extremely wide measuring range and excellent resolution
- One and the same force plate covers broad range of applications
- High overload protection and virtually no deterioration with age

Kistler guarantees reliable specifications and calibration

- High accuracy, linearity and sensitivity
- High natural frequency and damping
- Worldwide sales and customer service network

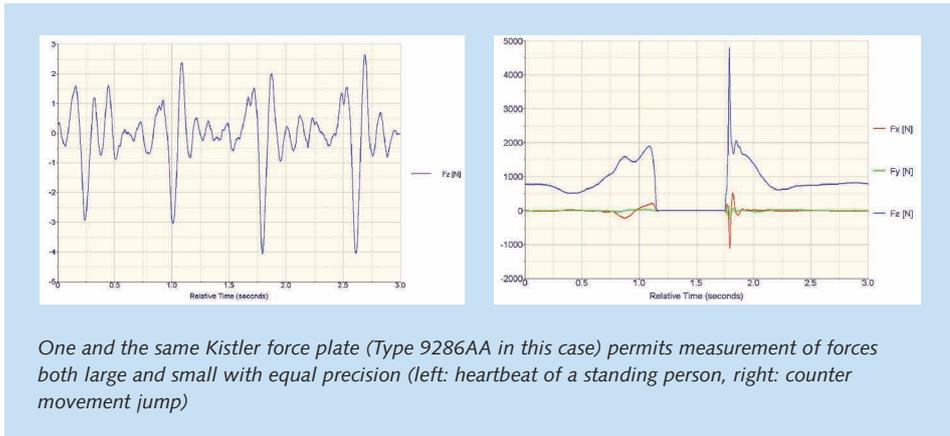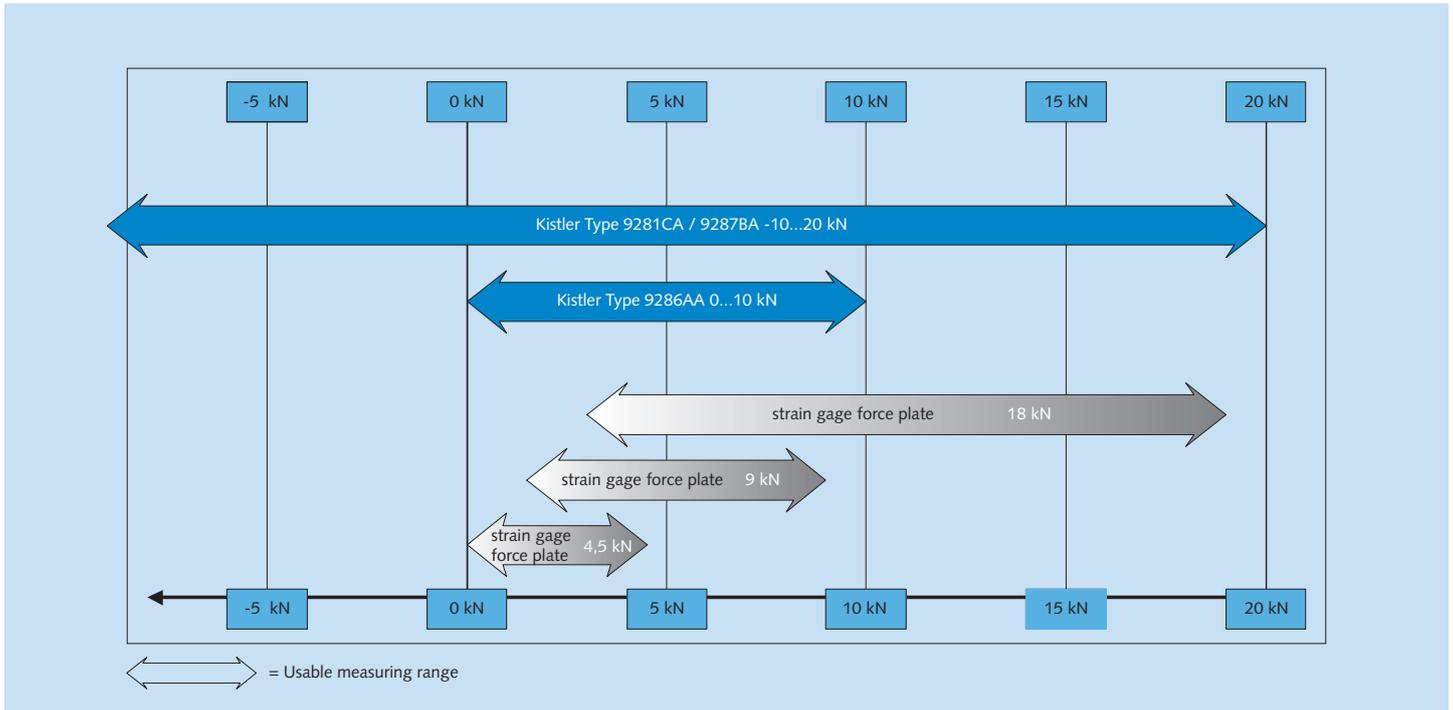

One and the same Kistler force plate (Type 9286AA in this case) permits measurement of forces both large and small with equal precision (left: heartbeat of a standing person, right: counter movement jump)

measure. analyze. innovate.

Kalibrierschein Kraft Calibration Certificate Force

Type: 9281CA Serial No. 1415655

Kalibriert durch Calibrator No. S. Magarinos	Datum 04.02.05
Referenzgenelle Referenzinstrument Gebrauchsinstrument (Modell / Seriennummer)	Type Typ Kistler 9062/08 Kistler 9062/08 Kistler 9062/08 Kistler 9281CA
Kontrollinstanz Control Unit	Series-No. Ser.-No.
Spannungskategorie Voltage, Calibration	138643 127874 474251 562527
Umgangstemperatur Ambient Temperature	SIEMENS 373656
Relative Feuchte Relative Humidity	
°C	%
°C	%

Skizze
Sketch

Messergebnisse Results of Measurement

Bereich Range	Kalibrierbereich Calibration Range		Überspannung Overrange							
	F_x / F_y	F_z	$F_x \rightarrow F_x$	$F_x \rightarrow F_y$	$F_y \rightarrow F_x$	$F_y \rightarrow F_y$	$F_x \rightarrow F_z$	$F_y \rightarrow F_z$	$F_z \rightarrow F_x$	$F_z \rightarrow F_y$
	mm	mm	%	%	%	%	%	%	%	%
1^{st}	0 - 5	0 - 20	0.02	0.02	0.1	-0.7	0.1	0.1	0.1	-0.2
2^{nd}	0 - 10	0 - 40	0.02	0.02	0.1	-0.7	0.1	0.1	0.1	-0.2
3^{rd}	0 - 20	0 - 80	0.02	0.02	0.1	-0.7	0.1	0.1	0.1	-0.2
4^{th}	0 - 50	0 - 200	0.02	0.02	0.1	-0.7	0.1	0.1	0.1	-0.2
5^{th}	0 - 100	0 - 400	0.02	0.02	0.1	-0.7	0.1	0.1	0.1	-0.2

1^{st} Greatest Absolute 2^{nd} Estimated Character No.
 3^{rd} Smallest Calculated (0.05%) 4^{th} Signal Signal

Messverfahren Konventionelle Kalibrierung, Vergleichsverfahren
Measurement Procedure Conventional Calibration, Comparison Method

Bestätigung Confirmation

Wir bestätigen, dass das oben identifizierte Gerät nach den vorgeschriebenen Verfahren geprüft wurde. Alle Messwerte sind auf nationale Normen rückverfolgbar. Kistler bietet die SCC-Classic Calibration Services Kalibrierstelle Nr. 049, akkreditiert nach ISO 17025. Das Kalibrationsmessungssystem ist nach ISO 9001 zertifiziert.

We confirm that the device identified above was tested by the prescribed procedures. All measuring devices are traceable to national standards. The SCC-Classic Calibration Services Calibration Laboratory No. 049 is operated by Kistler and accredited per ISO 17025. The calibration measurement system is certified per ISO 9001.

Kistler Instrumente AG Kistler Messtechnik GmbH Kistler Messtechnik GmbH Kistler Messtechnik GmbH
Säulstrasse 22 Tel. +41 52 221 11 11 D 78193 Wetzlar SC 732 Bank: CH17 0020 0110 3003 7402 8
Postfach Fax +41 52 221 14 14 Tele: 0049 042704 Mobile: No. 225 710
CH-8408 Wetzlar info@kistler.com Konto: 1122 0794 428 SO 9001 zertifiziert www.kistler.com
5610 exp. 1/1

Kistler worldwide

Europe

Austria

Kistler GmbH
Lemböckgasse 49f
1230 Wien
Tel. +43 1 867 48 67 0
sales.at@kistler.com

Czech Republic/Slovakia

Kistler, s.r.o.
Zelený pruh 99/1560
140 00 Praha 4
Tel. +420 296 374 878
sales.cz@kistler.com

Denmark/Norway/Sweden

Kistler Nordic AB
Aminogatan 34
431 53 Mölndal
Tel. +46 31 871 566
info.se@kistler.com

Finland

Kistler Nordic AB
Särkiniementie 3
00210 Helsinki
Tel. +358 9 612 15 66
sales.fi@kistler.com

France

Kistler France
ZA de Courtabœuf 1
15, avenue du Hoggar
91953 Les Ulis cedex
Tel. +33 1 69 18 81 81
info.fr@kistler.com

Germany

Kistler Instrumente GmbH
Daimlerstrasse 6
73760 Ostfildern
Tel. +49 711 34 07 0
info.de@kistler.com

Italy

Kistler Italia s.r.l.
Via Ruggero di Lauria, 12/B
20149 Milano
Tel. +39 02 481 27 51
sales.it@kistler.com

Netherlands

Kistler B.V. Nederland
Leeghwaterstraat 25
2811 DT Reeuwijk
Tel. +31 182 304 444
sales.nl@kistler.com

Switzerland/Liechtenstein

Kistler Instrumente AG
Verkauf Schweiz
Eulachstrasse 22
8408 Winterthur
Tel. +41 52 224 12 32
sales.ch@kistler.com

United Kingdom

Kistler Instruments Ltd.
13 Murrell Green Business Park
London Road
Hook, Hampshire RG27 9GR
Tel. +44 1256 74 15 50
sales.uk@kistler.com

Asia

Japan

Kistler Japan Co., Ltd.
23rd floor, New Pier Takeshiba North Tower
1-11-1, Kaigan, Minato-ku
Tokyo 105-0022
Tel. +81 3 3578 0271
sales.jp@kistler.com

China, People's Republic of

Kistler China Ltd.
Room 925, Yuan Chen Xin Building
No. 12 E1, Yuminlu Road Deshengmenwai
Beijing 100029
Tel. +86 10 8225 2163
sales.cn@kistler.com

India

Kistler Instruments (Pte) Ltd.
India Liaison Office
2B Century Plaza
560/562 Anna Salai
Teynampet, Chennai 600 018
Tel. +91 44 4213 2089
sales.in@kistler.com

Korea, Republic of

Kistler Korea Co., Ltd.
Gyeonggi Venture Anyang
Technical College Center 801
572-5, Anyang-Dong, Manan-Gu,
Anyang-City, Gyeonggi-Do 430-731
Tel. +82 31 465 6013
sales.kr@kistler.com

Singapore

Kistler Instruments (Pte) Ltd.
50 Bukit Batok Street 23
#04-06 Midview Building
Singapore 659578
Tel. +65 6316 7331
sales.sg@kistler.com

Taiwan

Kistler Representative Office in Taiwan
Room 9, 8F, No. 6, Lane 180
Sec. 6, Mincyuan E. Road
Taipei 114
Tel. +886 2 7721 2121
sales.tw@kistler.com

Thailand

Kistler Instrument (Thailand) Co., Ltd.
26/56 TPI Tower, 20th Floor
Nanglingee Rd., (Chan Tat Mai Rd.)
Thungmahamek, Sathorn
Bangkok 10120
Tel. +66 2678 6779-80
sales.thai@kistler.com

America

USA/Canada/Mexico

Kistler Instrument Corp.
75 John Glenn Drive
Amherst, NY 14228-2171
Tel. +1 716 691 5100
sales.us@kistler.com

Australia

Australia

Kistler Instruments Australia Pty Ltd
5 Glenn Court
Rowville, Victoria 3178
Tel. +61 3 9755 8222
sales.au@kistler.com

Other countries

Kistler Instrumente AG

Export Sales
Eulachstrasse 22, 8408 Winterthur
Switzerland
Tel. +41 52 224 11 11
sales.export@kistler.com

www.kistler.com

KISTLER
measure. analyze. innovate.

Headquarters

Switzerland

Kistler Group
Eulachstrasse 22, 8408 Winterthur
Tel. +41 52 224 11 11
Fax +41 52 224 14 14
info@kistler.com